VetLine Canine Distemper Virus (CDV)

ELISA

Enzyme immunoassay for the qualitative determination of antibodies against Canine Distemper Virus (CDV) in veterinary serum.

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ENGLISH

1. INTRODUCTION

Canine distemper is a contagious, febrile disease of canids and other carnivores caused by Canine Distemper Virus (CDV). CDV is a member of the Paramyxoviridae family and the morbillivirus genus.

Canine Distemper Virus is shed in all excretions and secretions, and is spread by direct and indirect contact via ingestion or inhalation of aerosolised droplets. Initial infection of the nasal cavity, pharynx and lungs is followed by a macrophage-associated viraemia to local lymph nodes where the virus replicates. Within one week, all lymphoid tissue is infected, causing lymphopenia. A second viraemic stage then distributes the virus to the surface epithelium of the respiratory, gastrointestinal and urogenital tracts, endocrine tissue and the grey and white matter of the CNS. A biphasic pyrexia is typical of distemper infection: the first fever occurs 3-6 days post-infection and is associated with lymphopenia, and the second peak coincides with widespread viraemia. Further signs depend on both the virus strain and the immune response mounted. In the event of a strong humoral and cellular response, disease may remain subclinical, and if a weak immune response is mounted infection is generally subacute. If the immune response fails, acute disease and potentially death ensues. When clinical disease manifests, this is initially characterized by lethargy, dehydration, anorexia, and weight loss followed by more specific signs depending on the principally affected organ.

Canine distemper is often fatal, but an increased production of virus-neutralizing antibodies can promote the recovery of the animal. However, CDV can persist in the uvea, CNS, lymphoid organs and footpads despite elimination from most organs and the blood. This can result in "old dog encephalitis" in dogs that recovered from acute canine distemper years previously. In this, several neurological episodes occur over weeks to months, and usually culminate in the death of the dog.

Pathogen	Disease	Symptomes
Canine Distemper Virus (CDV)	Canine distemper	<ul style="list-style-type: none">• Rhinitis• Conjunctivitis• Leukopenia• Pneumonia cough with purulent, bloody sputum• Tonsillitis• Enteritis• Secondary bacterial infections (for example, Bordetella)• Neurological symptoms such as tremor, ataxia, nystagmus, optic nerve damage• Pustular Dermatitis• Hyperkeratosis of the footpads and the nose leather (hard pad disease)• Enamel hypoplasia ("distemper teeth") in puppies

Infections may be diagnosed by:

- PCR
- Serology: Detection of antibodies by ELISA

2. INTENDED USE

The NovaTec VetLine Canine Distemper Virus (CDV) ELISA is intended for the qualitative determination of antibodies against Canine Distemper Virus (CDV) in veterinary serum.

3. PRINCIPLE OF THE ASSAY

The qualitative immunoenzymatic determination of specific antibodies is based on the ELISA (Enzyme-linked Immunosorbent Assay) technique.

Microtiterplates are coated with specific antigens to bind corresponding antibodies of the sample. After washing the wells to remove all unbound sample material a horseradish peroxidase (HRP) labelled conjugate is added. This conjugate binds to the captured antibodies. In a second washing step unbound conjugate is removed. The immune complex formed by the bound conjugate is visualized by adding Tetramethylbenzidine (TMB) substrate which gives a blue reaction product.

The intensity of this product is proportional to the amount of specific antibodies in the sample. Sulphuric acid is added to stop the reaction. This produces a yellow endpoint colour. Absorbance at 450/620 nm is read using an ELISA Microtiterplate reader.

4. MATERIALS

4.1. Reagents supplied

- **Microtiterplate:** 12 break-apart 8-well snap-off strips coated with Canine Distemper Virus (CDV) antigens; in resealable aluminium foil.
- **Sample Dilution Buffer:** 1 bottle containing 100 mL of phosphate buffer (10 mM) for sample dilution; pH 7.2 ± 0.2; coloured yellow; ready to use; white cap; ≤ 0.0015 % (v/v) CMIT/MIT (3:1).
- **Stop Solution:** 1 bottle containing 15 mL sulphuric acid, 0.2 mol/L; ready to use; red cap.
- **Washing Buffer (20x conc.):** 1 bottle containing 50 mL of a 20-fold concentrated phosphate buffer (0.2 M), pH 7.2 ± 0.2, for washing the wells; white cap.
- **Conjugate:** 1 bottle containing 20 mL of peroxidase labelled Protein A/G; coloured yellow; ready to use; white cap; ≤ 0.02 % (v/v) MIT.
- **TMB Substrate Solution:** 1 bottle containing 15 mL 3,3',5,5'-tetramethylbenzidine (TMB), < 0.1 %; ready to use; yellow cap.
- **Positive Control:** 1 vial containing 2 mL; coloured yellow; ready to use; red cap; ≤ 0.02 % (v/v) MIT.
- **Cut-off Control:** 1 vial containing 3 mL; coloured yellow; ready to use; green cap; ≤ 0.02 % (v/v) MIT.
- **Negative Control:** 1 vial containing 2 mL; coloured yellow; ready to use; blue cap; ≤ 0.0015 % (v/v) CMIT/ MIT (3:1).

For hazard and precautionary statements see 12.1

For potential hazardous substances please check the safety data sheet.

4.2. Materials supplied

- 1 Cover foil
- 1 Instruction for use (IFU)
- 1 Plate layout

4.3. Materials and Equipment needed

- ELISA Microtiterplate reader, equipped for the measurement of absorbance at 450/620 nm
- Incubator 37 °C
- Manual or automatic equipment for rinsing Microtiterplate wells
- Pipettes to deliver volumes between 10 and 1000 µL
- Vortex tube mixer
- Distilled water
- Disposable tubes

5. STABILITY AND STORAGE

Store the kit at 2...8 °C. The opened reagents are stable up to the expiry date stated on the label when stored at 2...8 °C.

6. REAGENT PREPARATION

It is very important to bring all reagents and samples to room temperature (20...25 °C) and mix them before starting the test run!

6.1. Microtiterplate

The break-apart snap-off strips are coated with Canine Distemper Virus (CDV) antigens. Immediately after removal of the strips, the remaining strips should be resealed in the aluminium foil along with the desiccant supplied and stored at 2...8 °C.

6.2. Washing Buffer (20x conc.)

Dilute Washing Buffer 1 + 19; e. g. 10 mL Washing Buffer + 190 mL distilled water. The diluted buffer is stable for 5 days at room temperature (20...25 °C). In case crystals appear in the concentrate, warm up the solution to 37 °C e.g. in a water bath. Mix well before dilution.

6.3. TMB Substrate Solution

The reagent is ready to use and has to be stored at 2...8 °C, away from the light. The solution should be colourless or could have a slight blue tinge. If the substrate turns into blue, it may have become contaminated and should be thrown away.

7. SAMPLE COLLECTION AND PREPARATION

Use canine serum samples with this assay. If the assay is performed within 5 days after sample collection, the samples should be kept at 2...8 °C; otherwise they should be aliquoted and stored deep-frozen (-70...-20 °C). If samples are stored frozen, mix thawed samples well before testing. Avoid repeated freezing and thawing. Heat inactivation of samples is not recommended.

7.1. Sample Dilution

Before assaying, all samples should be diluted 1+100 with Sample Dilution Buffer. Dispense 10 µL sample and 1 mL Sample Dilution Buffer into tubes to obtain a 1+100 dilution and thoroughly mix with a Vortex.

8. ASSAY PROCEDURE

Please read the instruction for use carefully **before** performing the assay. Result reliability depends on strict adherence to the instruction for use as described. The following test procedure is only validated for manual procedure. If performing the test on ELISA automatic systems we recommend increasing the washing steps from three up to five and the volume of Washing Buffer from 300 μ L to 350 μ L to avoid washing effects. Pay attention to chapter 12. Prior to commencing the assay, the distribution and identification plan for all samples and standards/controls (duplicates recommended) should be carefully established on the plate layout supplied in the kit. Select the required number of microtiter strips or wells and insert them into the holder.

Perform all assay steps in the order given and without any delays.

A clean, disposable tip should be used for dispensing each standard/control and sample.

Adjust the incubator to 37 ± 1 °C.

1. Dispense 100 μ L standards/controls and diluted samples into their respective wells. Leave well A1 for the Substrate Blank.
2. Cover wells with the foil supplied in the kit.
3. **Incubate for 1 hour \pm 5 min at 37 ± 1 °C.**
4. When incubation has been completed, remove the foil, aspirate the content of the wells and wash each well three times with 300 μ L of Washing Buffer. Avoid overflows from the reaction wells. The interval between washing and aspiration should be > 5 sec. At the end carefully remove remaining fluid by tapping strips on tissue paper prior to the next step!
Note: Washing is important! Insufficient washing results in poor precision and false results.
5. Dispense 100 μ L Conjugate into all wells except for the Substrate Blank well A1.
6. **Incubate for 30 min at room temperature (20...25 °C).** Do not expose to direct sunlight.
7. Repeat step 4.
8. Dispense 100 μ L TMB Substrate Solution into all wells.
9. **Incubate for exactly 15 min at room temperature (20...25 °C) in the dark.** A blue colour occurs due to an enzymatic reaction.
10. Dispense 100 μ L Stop Solution into all wells in the same order and at the same rate as for the TMB Substrate Solution, thereby a colour change from blue to yellow occurs.
11. Measure the absorbance at 450/620 nm within 30 min after addition of the Stop Solution.

8.1. Measurement

Adjust the ELISA Microtiterplate reader **to zero** using the **Substrate Blank**.

If - due to technical reasons - the ELISA Microtiterplate reader cannot be adjusted to zero using the Substrate Blank, subtract its absorbance value from all other absorbance values measured in order to obtain reliable results!

Measure the absorbance of all wells at **450 nm** and record the absorbance values for each standard/control and sample in the plate layout.

Bichromatic measurement using a reference wavelength of 620 nm is recommended.

Where applicable calculate the **mean absorbance values** of all duplicates.

9. RESULTS

9.1. Run Validation Criteria

In order for an assay to be considered valid, the following criteria must be met:

- **Substrate Blank:** Absorbance value < 0.100
- **Negative Control:** Absorbance value < 0.200 and $< \text{Cut-off}$
- **Cut-off Control:** Absorbance value $0.150 - 1.300$
- **Positive Control:** Absorbance value $> \text{Cut-off}$

If these criteria are not met, the test is not valid and must be repeated.

9.2. Calculation of Results

The Cut-off is the mean absorbance value of the Cut-off Control determinations.

Example: Absorbance value Cut-off Control 0.44 + absorbance value Cut-off control 0.42 = 0.86 / 2 = 0.43
Cut-off = 0.43

9.2.1. Results in Units [NTU]

$$\frac{\text{Sample (mean) absorbance value} \times 10}{\text{Cut-off}} = [\text{NovaTec Units} = \text{NTU}]$$

Example:
$$\frac{1.591 \times 10}{0.43} = 37 \text{ NTU}$$

9.3. Interpretation of Results

Normal value ranges for this ELISA should be established by each laboratory based on its own sample populations in the geographical areas serviced.

The following values should be considered as a guideline:

Cut-off	6,5 NTU	-
Positive	> 7 NTU	Antibodies against the pathogen are present. There has been a contact with the antigen (pathogen resp. vaccine).
Equivocal	6 – 7 NTU	Antibodies against the pathogen could not be detected clearly. It is recommended to repeat the test with a fresh sample in 2 to 4 weeks. If the result is equivocal again the sample is judged as negative .
Negative	< 6 NTU	The sample contains no antibodies against the pathogen. A previous contact with the antigen (pathogen resp. vaccine) is unlikely.
Diagnosis of an infectious disease should not be established on the basis of a single test result. A precise diagnosis should take into consideration clinical history, symptomatology as well as serological data.		

10. SPECIFIC PERFORMANCE CHARACTERISTICS

The results refer to the groups of samples investigated; these are not guaranteed specifications.

The performance data have been established with selected serum samples. Due to the nature of the Protein A/G conjugate this ELISA should react with other mammalian species also. More detailed information is available on request.

10.1. Precision

Intraassay	n	Mean (E)	CV (%)
#1	24	0.314	6.93
#2	24	0.931	3.70
#3	24	0.718	9.98

Interassay	n	Mean (NTU)	CV (%)
#1	12	33.93	2.76
#2	12	28.09	7.46
#3	12	3.85	12.11

10.2. Diagnostic Specificity

The diagnostic specificity is defined as the probability of the assay of scoring negative in the absence of the specific analyte.
Diagnostic Specificity canine: 92.31 % (95 % confidence interval: 63.97 % - 99.81 %)

10.3. Diagnostic Sensitivity

The diagnostic sensitivity is defined as the probability of the assay of scoring positive in the presence of the specific analyte.
Diagnostic Sensitivity canine: 96.88 % (95 % confidence interval: 83.78 % - 99.92 %)

10.4. Interferences

Interferences with hemolytic, lipemic or icteric samples are not observed up to a concentration of 10 mg/mL hemoglobin, 5 mg/mL triglycerides and 0.5 mg/mL bilirubin.

10.5. Cross Reactivity

Cross reactions cannot be excluded.

11. LIMITATIONS OF THE PROCEDURE

Bacterial contamination or repeated freeze-thaw cycles of the sample may affect the absorbance values.

12. PRECAUTIONS AND WARNINGS

- Only for veterinary in-vitro diagnostic use.
- All materials of human or animal origin should be regarded and handled as potentially infectious.
- All components of human origin used for the production of these reagents have been tested for anti-HIV antibodies, anti-HCV antibodies and HBsAg and have been found to be non-reactive.
- Do not interchange reagents or strips of different production lots.
- No reagents of other manufacturers should be used along with reagents of this test kit.
- Do not use reagents after expiry date stated on the label.
- Use only clean pipette tips, dispensers, and lab ware.
- Do not interchange screw caps of reagent vials to avoid cross-contamination.
- Close reagent vials tightly immediately after use to avoid evaporation and microbial contamination.
- After first opening and subsequent storage check conjugate and standard/control vials for microbial contamination prior to further use.
- To avoid cross-contamination and falsely elevated results pipette patient samples and dispense reagents without splashing accurately into the wells.
- The ELISA is only designed for qualified personnel who are familiar with good laboratory practice.

12.1. Safety note for reagents containing hazardous substances

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Therefore, the following hazard and precautionary statements apply.

Warning

H317	May cause an allergic skin reaction.
P261	Avoid breathing spray
P280	Wear protective gloves/protective clothing.
P302+P352	IF ON SKIN: Wash with plenty of soap and water.
P333+P313	If skin irritation or rash occurs: Get medical advice/attention.
P362+P364	Take off contaminated and Wash it before reuse.

Further information can be found in the safety data sheet.

12.2. Disposal Considerations

Residues of chemicals and preparations are generally considered as hazardous waste. The disposal of this kind of waste is regulated through national and regional laws and regulations. Contact your local authorities or waste management companies which will give advice on how to dispose hazardous waste.

13. ORDERING INFORMATION

Prod. No.: MEAVT0330 VetLine Canine Distemper Virus (CDV) ELISA (96 Determinations)

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Tilley, L.P. and Smith, F.W.K.(2004)The 5-minute Veterinary Consult (Fourth Edition) Blackwell Publishing.

Merck & Co (2008) The Merck Veterinary Manual (Eighth Edition) Merial.

Carter, GR and Wise, DJ (2005) A Concise Review of Veterinary Virology, International Veterinary Information Service.

Vandevelde, M and Zurbriggen, A (1995) The neurobiology of canine distemper virus infection. Veterinary Microbiology 44, 271-280.

Chappuis, G (1995) Control of canine distemper. Veterinary Microbiology, 44, 351-358.

Beineke, A, Puff, C, Seehusen, F, Baumgartner, W (2009) Pathogenesis and immunopathology of systemic and nervous canine distemper. Veterinary Immunology and Immunopathology, 127, 1-18.

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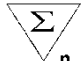
SCHEME OF THE ASSAY

VetLine Canine Distemper Virus (CDV) ELISA

Test Preparation

Prepare reagents and samples as described.
 Establish the distribution and identification plan for all samples and standards/controls on the plate layout supplied in the kit.
 Select the required number of microtiter strips or wells and insert them into the holder.

Assay Procedure

	Substrate Blank (A1)	Negative Control	Cut-off Control	Positive Control	Sample (diluted 1+100)
Negative Control	-	100 µL	-	-	-
Cut-off Control	-	-	100 µL	-	-
Positive Control	-	-	-	100 µL	-
Sample (diluted 1+100)	-	-	-	-	100 µL
Cover wells with foil supplied in the kit Incubate for 1 h at 37± °C Wash each well three times with 300 µL of Washing Buffer					
Conjugate	-	100 µL	100 µL	100 µL	100 µL
Incubate for 30 min at room temperature (20...25 °C) Do not expose to direct sunlight Wash each well three times with 300 µL of Washing Buffer					
TMB Substrate Solution	100 µL	100 µL	100 µL	100 µL	100 µL
Incubate for exactly 15 min at room temperature (20...25 °C) in the dark					
Stop Solution	100 µL	100 µL	100 µL	100 µL	100 µL
Photometric measurement at 450 nm (reference wavelength: 620 nm)					

NovaTec Immundiagnostica GmbH

Waldstraße 23 A6
 63128 Dietzenbach, Germany

Tel.: +49 (0) 6074-48760
 Email: info@NovaTec-ID.com
 Internet: www.NovaTec-ID.com

Fax: +49 (0) 6074-487629